



CMAXS

Installation/Un-installation

Manual






Be sure to login as administrator authority into your PC.

<INDEX>

1. Installation Manual	---2
2. Un-installation Manual	---5

1. Installation Manual


(1) Set CMAXS installation CD into your PC.

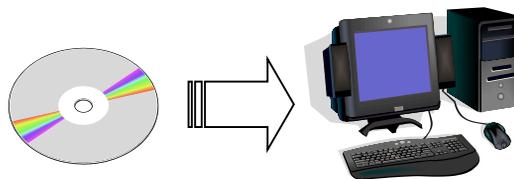

(2) "AutoPlay" screen is displayed.
Select [Run Launcher_PMS.exe].
(If it is not displayed, double click "Launcher_PMS.exe" in CD drive.)




(3) "CMAXS Setup Launcher" screen is automatically displayed.
Select [Install Program].

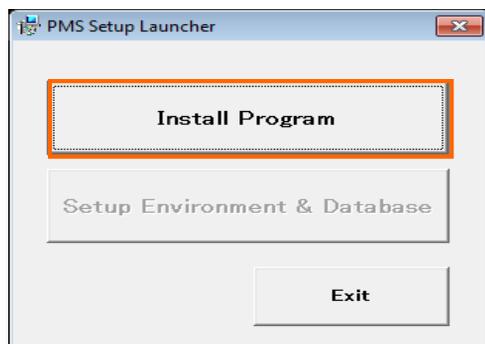


(4) "Microsoft .NET Framework 4" screen is displayed
Select [Accept].
The installation will start.

If this screen is not displayed, please skip to step (7)

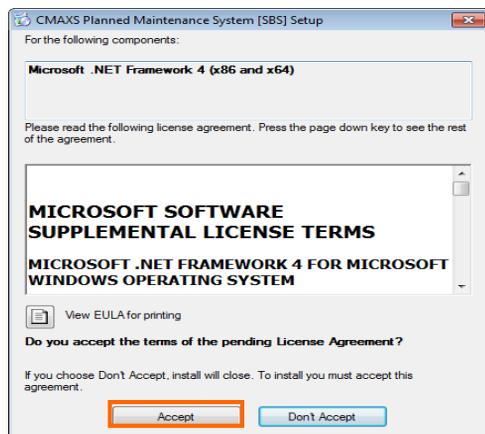


(5) After the installation is complete,
Restart screen is displayed.
Select [Yes].

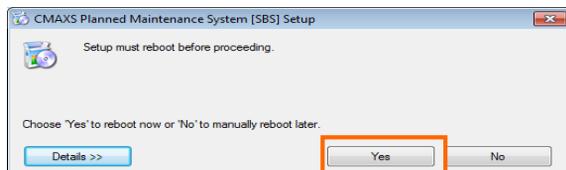


(6) After restart, “Setup Wizard” screen is displayed. Select [Next].




(7) “Select Installation Folder” screen is displayed. Make sure that check mark is put on “Everyone”, and then, select [Next].

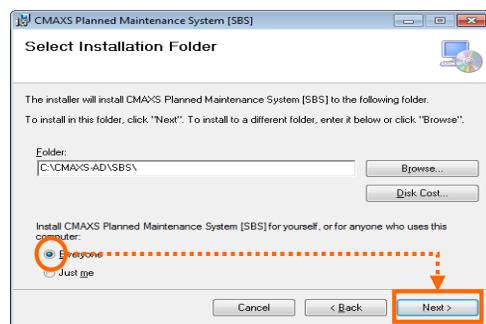

(8) “Confirm Installation” screen is displayed. Select [Next]. “Installing” screen is displayed, will start the installation.

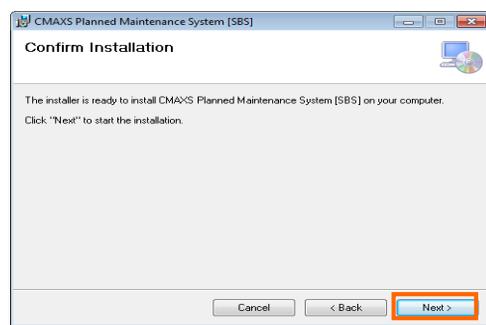


(9) After the installation is complete, “Installation Complete” screen is displayed. Select [Close].

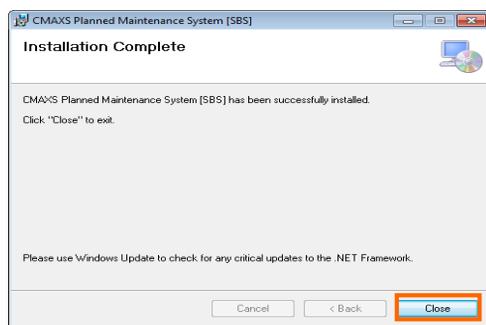


Please go to step (12).
If step (12) screen is not displayed, please go to step (10).


(10) Select start menu button, select [Computer].

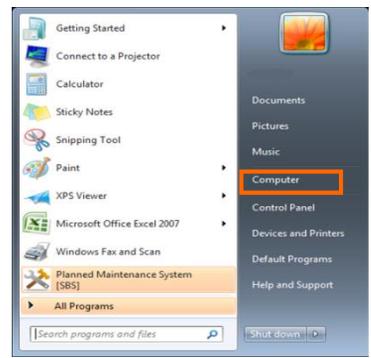

(11) Double click [DVD/CD-RW].

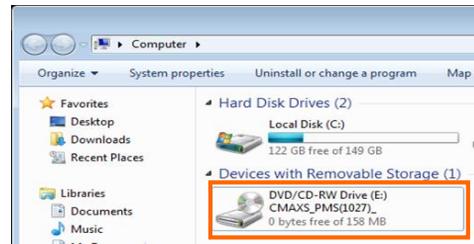


(12) "CMAXS Setup Launcher" screen is displayed.
Select [Setup Environment & Database].

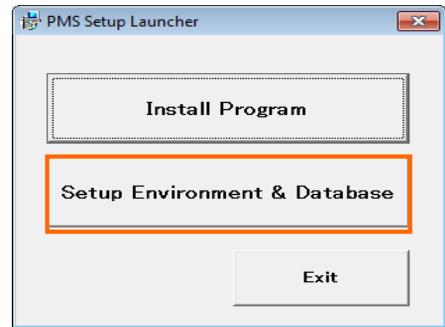


(13) "Setup Database & Environment" screen is displayed.
Make sure that check mark is put on [Installation of new Database],
and then, select [OK].




(14) "Is it correct in this setting" dialog is displayed.
Select [Yes].

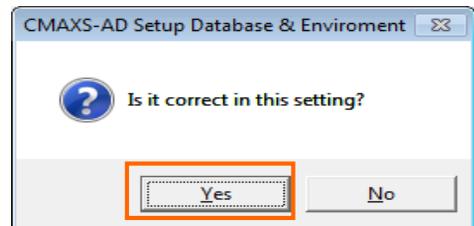

(15) "Database setup has been successfully completed" dialog is displayed.
Select [OK].

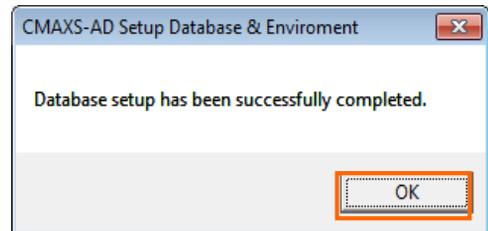

(16) "CMAXS Setup Launcher" screen is displayed.
Select [Exit].

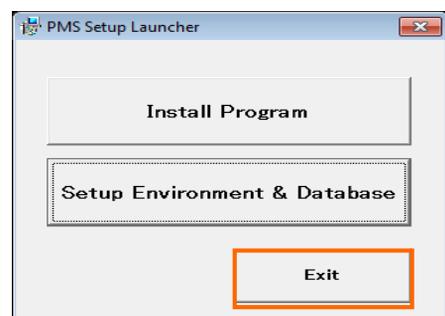


(17) CMAXS installation has been completed.
Please execute installation of the other system.

2. Un-installation Manual


(1) Select start menu bottom, select [Control Panel].




(2) Select [Program and features].

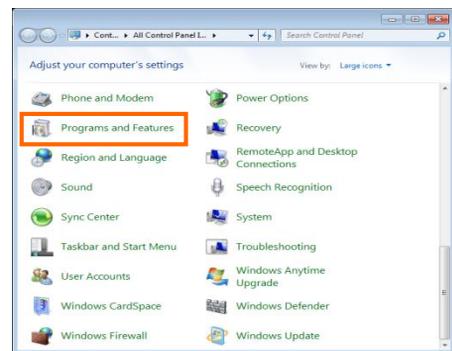

(3) "Uninstall or change a program" screen is displayed. Select [CMAXS system], and then, select [Uninstall].

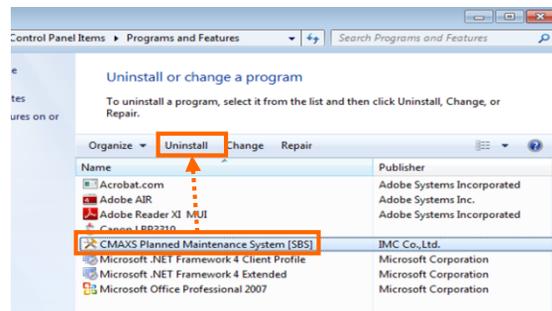

(4) Uninstall confirmation screen is displayed. Select [Yes].

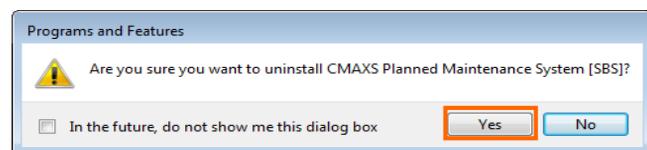

(5) Uninstallation is complete.